


## Capsule Filling Machine

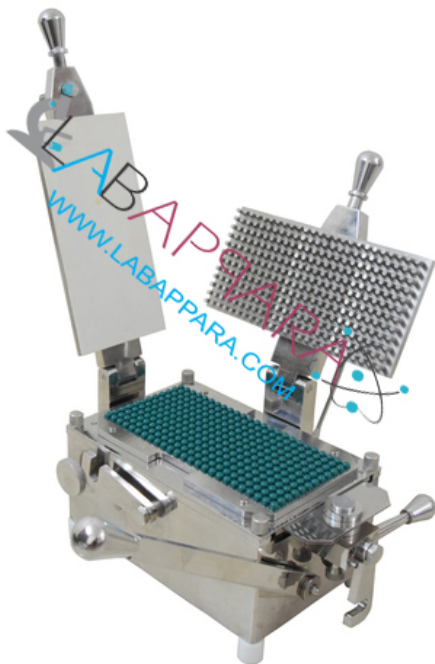

**Product Categories:** [Laboratory Equipments](#), [Pharmaceutical Equipments](#)

**Product Tags:** [biology lab equipment](#), [chemistry lab instruments](#)

**Product Page:** <https://www.labappara.com/product/capsule-filling-machine/>

### Product Description

#### Capsule Filling Machine - Scientific Equipment

This capsule filling performs a number of tasks that increases its demand in pharmaceutical industry. Semi Automatic Capsule Filling Machine consists of no. of plates for putting capsules by hand where as it also has function for removing caps of the capsules. Filling of powder in capsule is done by dispensing of brush pressing and after that fixing of caps is done.

A compact hand operated machine suitable to fill two-part hard gelatin capsules no. **0, 1 & 2** by interchangeable plates.

Main body made of heavy cast iron and all contact parts made of brass + richly plated or stainless steel.

Handle locking arrangement makes the operator to lift the loading plate without

fear of damaging the capsules.

Pressing is achieved by cam system which also saves time to press the plunger as compared to screw system.

**Features:**

Powder coated M.S body and all leavers

Stainless steel compressing pins plates for filling the capsules

All capsules touching parts in stainless steel

Easily lock the capsules remove the cap

Filling tray and brush is provided with machine

Heavy duty and manually operated

**Features:**

Energy efficient

Excellent functionality

Low operational noise