


CONSTANT LEVEL TANK

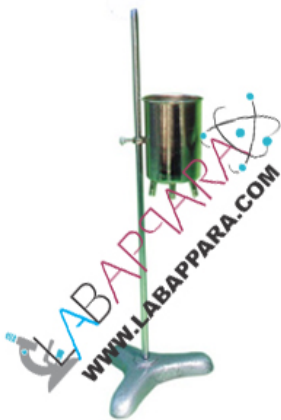

Product Categories: [Heat And Thermodynamics](#), [Physics Equipments](#)

Product Tags: [physics lab instruments](#)

Product Page: <https://www.labappara.com/product/constant-level-tank/>

Product Description

Constant Level Tank - Physics Equipment

For maintaining a uniform flow of water through calorimeter or other equipment

Regardless of pressure fluctuations in water supply line.

The level cannot rise above a central overflow tube in the chamber.

Consisting of a copper vessel of 100×75 mm diameter with 3 tubes for water inlet, outlet and overflow with boss head and stand.

A clamp for attaching to any 9 mm support rod is provided.

Supplied with suitable retort stand.