


## Dropping Funnels

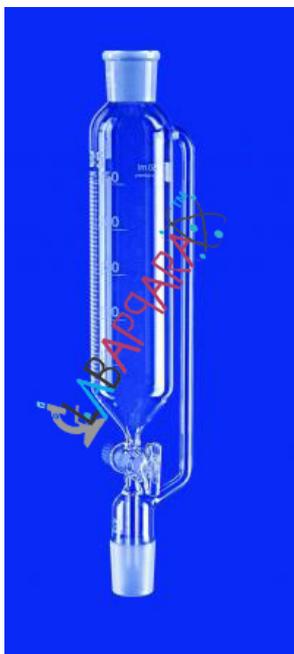

**Product Categories:** [Funnels Borosilicate Glass](#), [Laboratory Glassware](#)

**Product Tags:** [glassware](#), [Laboratory Equipment](#)

**Product Page:** <https://www.labappara.com/product/dropping-funnels/>

### Product Description

#### Dropping Funnels - Laboratory Glassware Equipment

"LABAPPARA", Cylindrical, with Interchangeable Stopper and Glass Key Stopcock, Plain Stem.

A dropping funnel is a type of laboratory glassware used to transfer fluids. They are fitted with a stopcock which allows the flow to be controlled. Dropping funnels are useful for adding reagents slowly, i.e. drop-wise. This may be desirable when the quick addition of the reagent may result in side reactions, or if the reaction is too vigorous.

Dropping funnels are usually fitted with a ground glass joint which allows the funnel to fit snugly onto, e.g. a round bottom flask. This also means it need not be clamped separately.

Capacity -50 ml, 100 ml, 150 ml, 250 ml, 500 ml, 1000 ml,