


## Laboratory Mill Willey Type

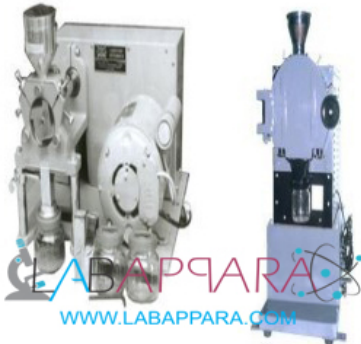

**Product Categories:** [Laboratory Equipments](#), [Pharmaceutical Equipments](#)

**Product Tags:** [biology lab equipment](#), [chemistry lab instruments](#)

**Product Page:** <https://www.labappara.com/product/laboratory-mill-willey-type/>

### Product Description

Laboratory Mill Willey Type - Scientific Equipment

For preparation of samples, for analysis with minimum loss of moisture of volatile constituents.

Weight distribution and balancing of rotating parts has been adjusted to provide greater stability and quieter performance with reduced vibration.

Unit having two stationary steel blades and Rotary with steel edges, provided in the cast Iron chamber.

A removable safety Glass sheet plate from the face of chamber of permits observation during operation.

Willey Mill is fitted on substantial Cast iron base.

Suitability for wider variety of commercial materials. fitted with 1/4 HP motor, supplied with nickel plated hopper of 50 mm diameter with cover, three delivery sieves

Having tops of 20, 40 and 80 mesh respectively, with plastic receiver, brush, L-N key, cord and plug.

Suitable to operate on 220 V, 1 Phase, 50 Cycles AC supply.

Complete unit supplied with aluminum / Cast chamber , polished hopper with cover,  
wooden plugged